



2006-10-3 Nanogold-Ab conjugates for fumonisin B1 analysis

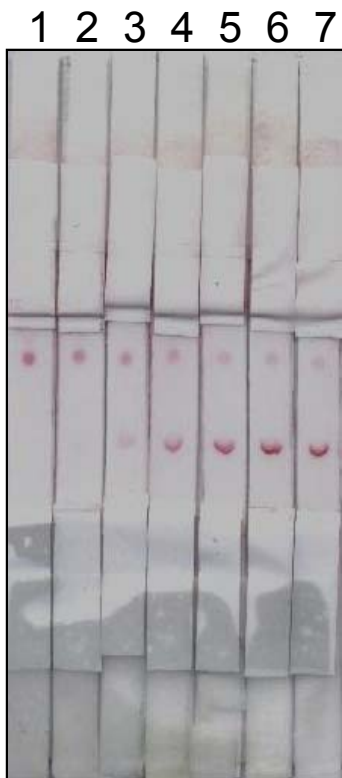

5 ul strip	1	2	3	4	5	6	7
Control line	Anti-R 1:10 1ul	Anti-R 1:10 1ul	Anti-R 1:10 1ul	Anti-R 1:10 1ul	Anti-R 1:10 1ul	Anti-R 1:10 1ul	Anti-R 1:10 1ul
Test line	FB1-OV 1ul	FB1-OV 1ul	FB1-OV 1ul	FB1-OV 1ul	FB1-OV 1ul	FB1-OV 1ul	FB1-OV 1ul
Gold conjugate Ab	FB1-40nm 5ul	FB1-40nm 5ul	FB1-40nm 5ul	FB1-40nm 5ul	FB1-40nm 5ul	FB1-40nm 5ul	FB1-40nm 5ul
sample	FB1 1ug/ml 200ul	FB1 100ng/ml 200ul	FB1 10ng/ml 200ul	FB1 1 ng/ml 200ul	FB1 0.1ng/ml 200ul	FB1 0.01ng/ml 200ul	0.01MPBS 200ul

2006-9-29 Nanogold-Ab conjugates for domoic acid analysis

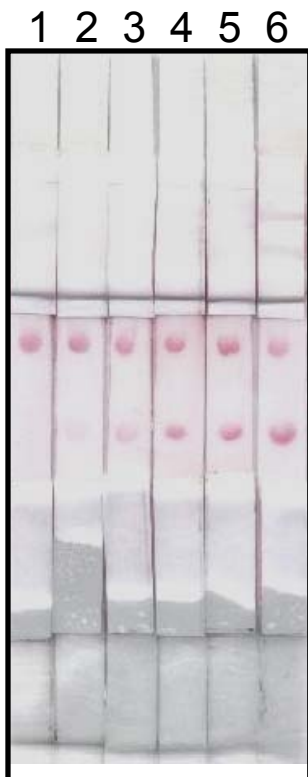

strip	1	2	3	4	5	6
Control line	Anti-R 1:1 0.5ul	Anti-R 1:1 0.5ul	Anti-R 1:1 0.5ul	Anti-R 1:1 0.5ul	Anti-R 1:1 0.5ul	Anti-R 1:1 0.5ul
Test line	DA-OVA 1ul	DA-OVA 1ul	DA-OVA 1ul	DA-OVA 1ul	DA-OVA 1ul	DA-OVA 1ul
Glod conjugate Ab	Yu4 1:10 5ul	Yu4 1:10 5ul	Yu4 1:10 5ul	Yu4 1:10 5ul	Yu4 1:10 5ul	Yu4 1:10 5ul
sample	DA 10ug/ml 200ul	DA 1ug/ml 200ul	DA 100ng/ml 200ul	DA 10ng/ml 200ul	DA 1ng/ml 200ul	0.01MPBS 200ul